

AD-A061 852

ARMY MATERIALS AND MECHANICS RESEARCH CENTER WATERTO--ETC F/G 13/11
BIAXIAL TESTING TECHNIQUES OF THIN-WALLED TUBULAR SPECIMENS.(U)
SEP 78 J J RAINY, R A SWANSON, S CHOU

UNCLASSIFIED


AMMRC-TR-78-37

NL

1 OF 1

AD
A061852



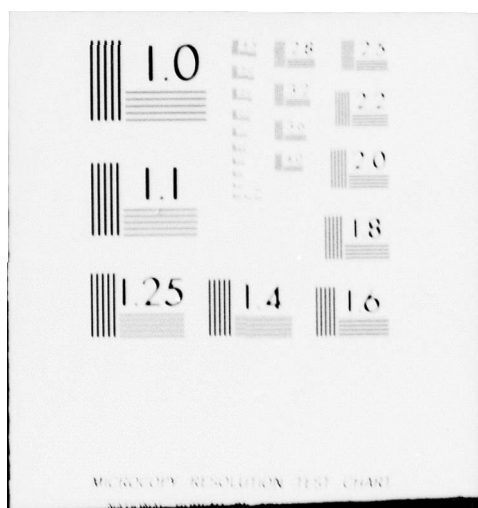

AMMRC TR 78-37

12 LEVEL II

AD

AD A061852

DDC FILE COPY

BIAXIAL TESTING TECHNIQUES OF THIN-WALLED TUBULAR SPECIMENS

JAMES H. RAINEY, RONALD A. SWANSON, and SHUN-CHIN CHOU
BALLISTIC MISSILE DEFENSE MATERIALS PROGRAM OFFICE

September 1978

Approved for public release; distribution unlimited.



ARMY MATERIALS AND MECHANICS RESEARCH CENTER
Watertown, Massachusetts 02172

78 12.1 016

The findings in this report are not to be construed as an official Department of the Army position, unless so designated by other authorized documents.

Mention of any trade names or manufacturers in this report shall not be construed as advertising nor as an official indorsement or approval of such products or companies by the United States Government.

DISPOSITION INSTRUCTIONS

Destroy this report when it is no longer needed.
Do not return it to the originator.

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

REPORT DOCUMENTATION PAGE		READ INSTRUCTIONS BEFORE COMPLETING FORM						
1. REPORT NUMBER AMMRC-TR-78-37	2. GOVT ACCESSION NO.	3. RECIPIENT'S CATALOG NUMBER						
4. TITLE (and Subtitle) BIAXIAL TESTING TECHNIQUES OF THIN-WALLED TUBULAR SPECIMENS	5. TYPE OF REPORT & PERIOD COVERED Final Report	6. PERFORMING ORG. REPORT NUMBER						
7. AUTHOR(s) James H. Rainey, Ronald A. Swanson Shun-Chin/Chou	8. CONTRACT OR GRANT NUMBER(s)							
9. PERFORMING ORGANIZATION NAME AND ADDRESS Army Materials and Mechanics Research Center Watertown, Massachusetts 02172 DRXMR-H	10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS D/A Project: SX3635040215 AMCMS Code: 633304.21500.03							
11. CONTROLLING OFFICE NAME AND ADDRESS U. S. Army Materiel Development and Readiness Command, Alexandria, Virginia 22333	12. REPORT DATE 11 Sep 1978	13. NUMBER OF PAGES 22						
14. MONITORING AGENCY NAME & ADDRESS (if different from Controlling Office) 12/28p	15. SECURITY CLASS. (of this report) Unclassified	15a. DECLASSIFICATION/DOWNGRADING SCHEDULE						
16. DISTRIBUTION STATEMENT (of this Report) Approved for public release; distribution unlimited.								
17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report)								
18. SUPPLEMENTARY NOTES								
19. KEY WORDS (Continue on reverse side if necessary and identify by block number) <table border="0"> <tr> <td>Biaxial stresses</td> <td>Fracture</td> </tr> <tr> <td>Mechanical tests</td> <td>Stress-strain</td> </tr> <tr> <td>Yield</td> <td>Aluminum</td> </tr> </table>			Biaxial stresses	Fracture	Mechanical tests	Stress-strain	Yield	Aluminum
Biaxial stresses	Fracture							
Mechanical tests	Stress-strain							
Yield	Aluminum							
20. ABSTRACT (Continue on reverse side if necessary and identify by block number) (SEE REVERSE SIDE) 403 105								

DDC
RECEIVED
DEC 6 1978
B

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE(When Data Entered)

Block No. 20

ABSTRACT

A biaxial testing technique of a thin-walled tubular specimen is described. The specimen was loaded with a combination of axial tension/compression and internal/external pressure. The technique can be used to test tubular specimens of any material. The yield surface of 2014-T651 aluminum was determined to illustrate the testing technique and data analysis.

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE(When Data Entered)

CONTENTS

	Page
I. INTRODUCTION.	1
II. TESTING TECHNIQUES.	1
A. Medium Strain Rate Machine.	1
B. Intensifier System.	2
C. Data Acquisition and Control System	2
D. Specimen Configuration and Grip Fixtures.	3
E. Specimen Strain Gaging.	6
F. Testing Control Parameters.	6
III. DATA ANALYSIS	6
IV. EXPERIMENTAL RESULTS AND DISCUSSIONS.	9
APPENDIX. BIAX AND BIANDV PROGRAMS	11

ACCESSION for	
NTIS	White Section <input checked="" type="checkbox"/>
DDC	Buff Section <input type="checkbox"/>
UNANNOUNCED	<input type="checkbox"/>
JUSTIFICATION	
BY	
DISTRIBUTION/AVAILABILITY CODES	
Dist. <input type="checkbox"/> and/or SPECIAL <input type="checkbox"/>	
A	

78 12:1 016

I. INTRODUCTION

The response of materials under multiaxial stress states beyond the elastic range is an essential part of the mechanical properties required by designers of structures. This report describes a technique to test material in the form of a thin-walled tube under axial and tangential (hoop) stresses. A brief description of the mechanical testing machine and information of an automated data acquisition and control system are presented. The automated control parameter determination and some experimental data of aluminum alloy 2014-T651 are also discussed.

II. TESTING TECHNIQUES

The equipment (Figure 1) required to perform the biaxial testing of tubular specimens can be described in three major components: namely, the medium strain rate machine, intensifier system, and data acquisition and control system.

The medium strain rate machine (MSRM) is used to generate the axial tensile or compression stress, while the intensifier system provides the tangential (or hoop) tensile or compression stress in a tubular specimen. The data acquisition and control system records and stores all data and generates command signals for both machines.

A. Medium Strain Rate Machine

The MSRM is a dual-mode test machine capable of generating 140,000 pounds axial tension or compression. The two modes of operation are an open loop system, and a closed loop servohydraulic system. The open loop system has the capability of strain rates from 1 to 50 sec^{-1} , but it is not normally used in biaxial testing. The closed loop system will produce strain rates from 10^{-5} to 10^{-1} sec^{-1} .

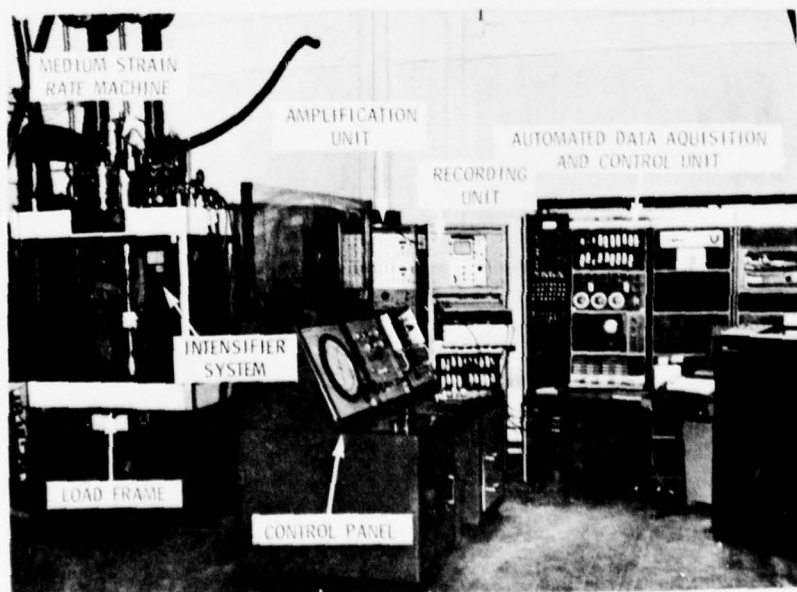


Figure 1. Automated materials characterization system.

The MSRM control panel gives the operator a selection of four different feedback control modes: load, displacement, strain, and optional. Therefore, with proper command inputs, tests at constant rates of load, displacement, and strain are performed. The operator may select one of many different load cells to achieve the best control over the desired testing range. The machine is equipped with a 15-gpm hydraulic power supply and a 15-gpm servovalve. Various fail-safe and limiting devices which either display a warning light or cause a machine shut-down are incorporated into the control system.

The MSRM load frame is designed for a stiffness greater than 15×10^6 lb/in. and has a total machine stretch of 0.005 inch at 140,000 lb load.

B. Intensifier System

The closed loop servohydraulic intensifier system used for the biaxial internal and external pressure testing is capable of generating a pressure of 100,000 psi. The intensifier vessel has a 1-inch inside diameter with an 8-inch length giving a volume of about 6 cubic inches. A Haskell air-driven pump is used to fill the intensifier system with Stoddard solvent as the pressurizing fluid. This fluid has a freezing point well in excess of seven kilobars.

The intensifier control panel gives the operator a selection of four different feedback loops: load, strain, displacement, and optional. With the proper command input, constant rates are possible under closed loop control. In essence, the intensifier system is a testing machine by itself. Different pressure transducers can be used so that the best control is available over the pressure range of interest. Pressure control is used in the optional mode. System pressure is also monitored using a precision Heise gage with a capacity of 100,000 psi. The intensifier uses the same hydraulic power supply as the medium strain rate machine.

C. Data Acquisition and Control System

The computer configuration consists of a central processor, 4K words of basic memory plus 12K words extended memory, real-time clock, relay register, 1.6 million word disk, two magnetic tape drives, teletype and line printer, display screen, multiplexed analog-to-digital converter (16 channels), and three digital-to-analog converters. The system interface includes eight active filters and scaling amplifiers.

The digital computer has a central processor which uses its memory to hold the operating system, to store programs during execution, and for temporary storage of data.

Command signals are generated by the central processor and sent to the digital-to-analog converters at a predetermined interval by the real-time clock. The converter changes binary numbers (12BITS), which are the internal information base of the computer, to an analog voltage ($\pm 10V$) that is acceptable to the servocontroller of the test machines. The controllers operate in a closed-loop mode generating a signal which drives the electrohydraulic servovalves. The servovalves regulate the flow of oil from the hydraulic power supply, to the test machine actuators, which deform the specimen. The specimen load, pressure, strains, and displacements

are continuously monitored and input into the computer. These analog signals are filtered to remove noise and scaled to match the input range of the analog-to-digital converter. The multiplexer selects one of the 16 channels for input to the analog-to-digital converter according to the program. A sample-and-hold circuit holds the value while it is converted to a binary number for input into the central processor. The data are stored in memory, then either displayed on the screen, printed, or stored on magnetic tape or disk. A second command signal may now be sent to other machines through a different digital-to-analog converter channel, based on the data received from the first signal. The real-time clock is used to determine sampling intervals, command intervals, and signal events. The teletype or console is used to input programs and parameters to the central processor. The magnetic tape units and disk provide fast access mass storage for programs and data. The relay register is used for additional control of the test machines and external equipment.

D. Specimen Configuration and Grip Fixtures

The thin-walled biaxial tubular specimen configuration shown in Figure 2 is designed for both internal and external pressure testing. The specimen designated as 601066-02 has a shorter gage length and overall length in comparison with the 601066-01 specimen. The shorter specimen is used to prevent premature buckling failure before the material reaches its yield stress when the specimen is subjected to compression in both axial and tangential directions.

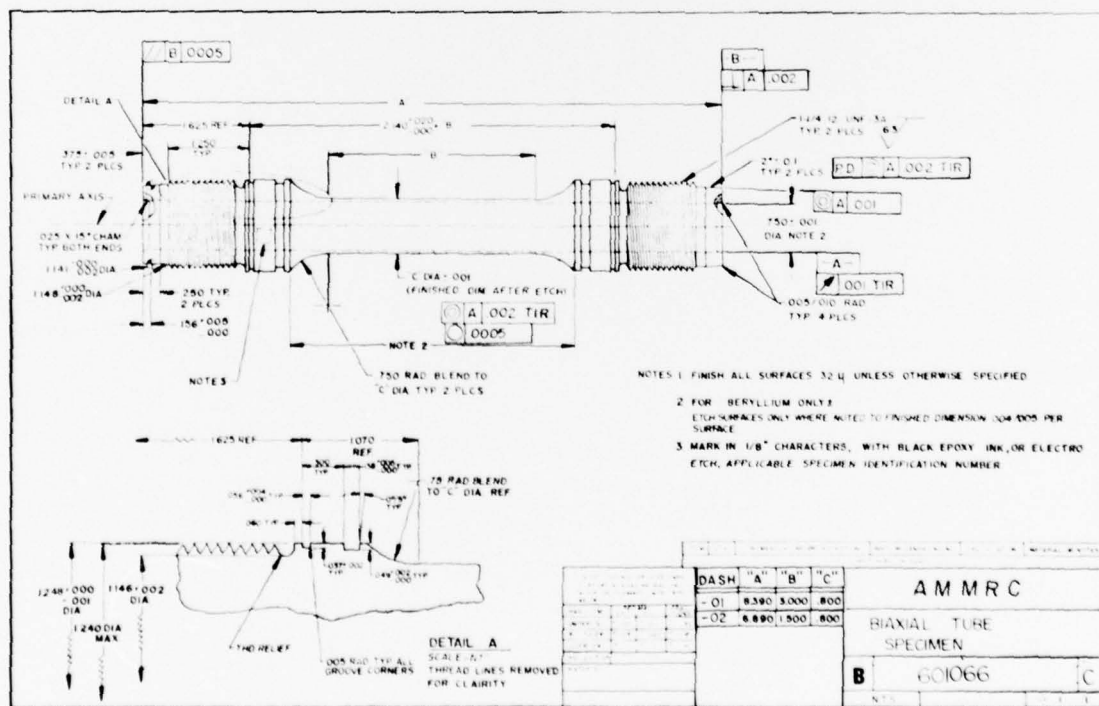


Figure 2. Specimen configuration for internal and external pressure tests.

In order to maintain the stress state in the specimen as close to a state of plane stress as possible, the wall thickness at the gage length is 0.025 inch, which results in an 8% difference in the tangential stress at the inner and outer radii.

The tensile tangential stress is created through the use of an internal plug assembly as shown in Figure 3, while the compression tangential stress is generated with a sleeve arrangement as shown in Figure 4. The external pressure assembly is shown in Figure 5. The pressurized fluid in the plug or sleeve is controlled by the intensifier system, and the pressure inside the plug or sleeve is contained by the use of "O" rings at both ends of the assembly.

The other precaution one must take in testing tubular specimens is that it is essential to maintain a homogeneous strain field in the gage section. This can be accomplished by a precise alignment procedure. A threaded collar containing eight bolt holes is threaded on each specimen end, having approximately 0.25 inch of the specimen exposed. The taper on the exposed ends of the specimen allows an initial alignment of the specimen with respect to the loading frame (see Figure 6). Bolts are then inserted through the collar and threaded to the load frame. While tightening bolts to the load frame, the two axial strain gages on the specimen are monitored to assure that, first, no bending moment is induced through uneven tightening of the bolts; and secondly, zero axial load is maintained. This procedure works very well as test results show that both axial strain gages record nearly identical strains until the specimen approaches failure.

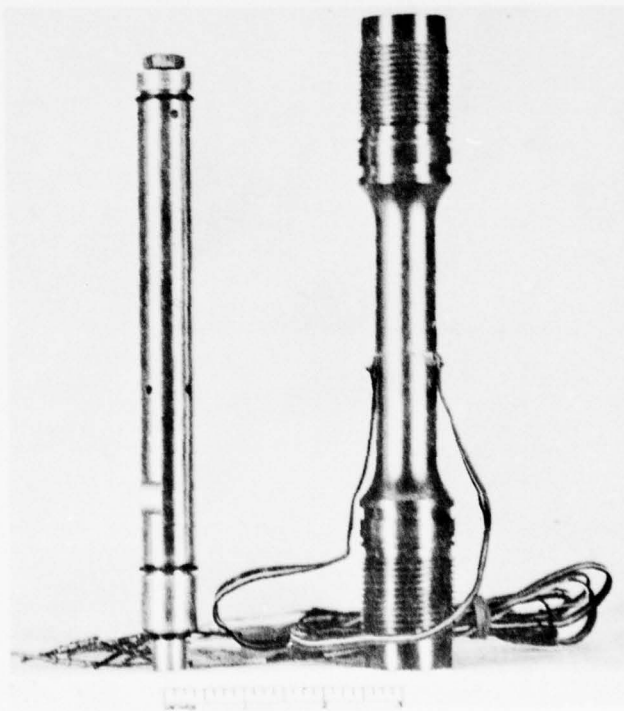


Figure 3. Biaxial internal plug and specimen.

19-066-7/AMC-78

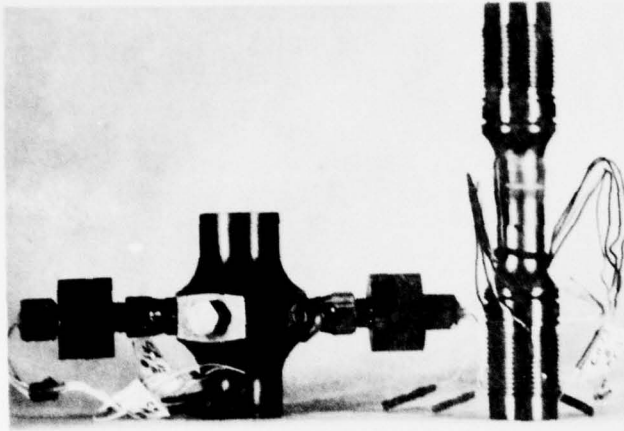


Figure 4. Biaxial external pressure sleeve and specimen.
19-066-6/AMC-78

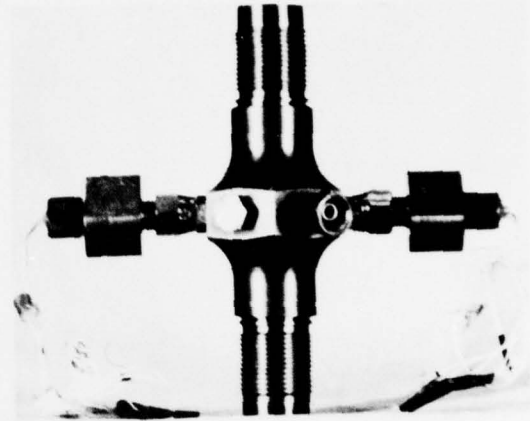


Figure 5. Biaxial external specimen assembly.
19-066-5/AMC-78

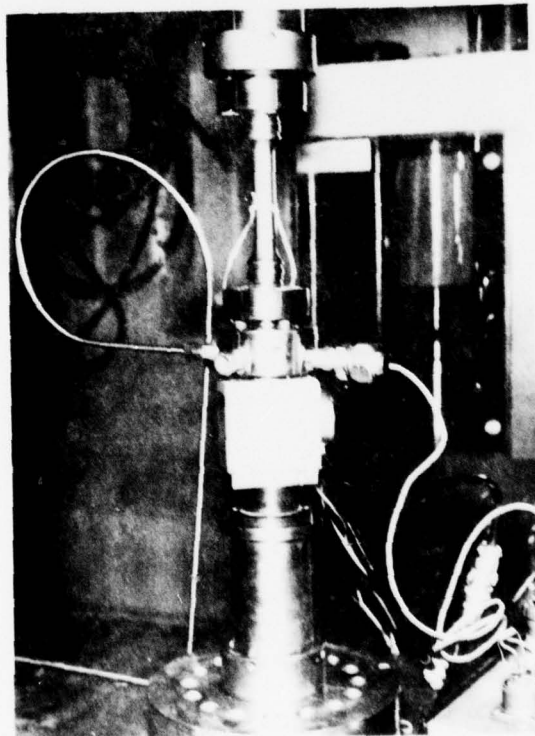


Figure 6. Grip arrangement for biaxial internal pressure test.

E. Specimen Strain Gaging

All biaxial specimens are instrumented with EP08-062TT-120, 90° "tee" rosette strain gages to measure the axial and tangential strains simultaneously. The gaging area is prepared in accordance with the manufacturer's instructions for the particular adhesive being used. Gages are bonded to the specimen with EPY-350 adhesive using a clamping pressure of 10 to 15 psi, and cured for two hours at 350 F. A protective coating of flowable, room temperature curing silicone rubber is used to protect the strain gages from environmental and handling damage. The detailed procedure for gaging specimens can be found in Micro Measurement Instruction Bulletin No. B-127-3, B-130-3, and B-137-2.

At the beginning of this study, four sets of "tee" rosette strain gages were used. They were located at the center of the gage section and placed 90° from each other along the circumference of the tubular specimen. However, since the alignment procedure described in the last paragraph creates a very good homogeneous strain field, it was decided that two sets of "tee" rosette strain gages are sufficient to measure the strain field in the tubular specimen. Test results reported are average values of the two sets of rosette strain gages located diametrically from each other.

F. Testing Control Parameters

Since the medium strain rate machine and intensifier system have separate servohydraulic mechanisms, the axial and tangential stresses (or strains) can, in principle, be controlled at different rates. This operation requires the use of a more sophisticated control theory because the material deformations are interrelated through the Poisson's effect. In this study, each test was carried out along a proportional load path which provides a constant effective strain rate in terms of either axial or tangential strain rate (see the mathematical formulation in the next section), that is, the ratio of the axial and tangential stresses is a constant throughout a test. As described in the section on "Data Acquisition and Control System," the signals from load cell and strain gages were magnified by the gage conditioning units; therefore, it would be logical to use the dominant strain as the control parameter. If the dominant strain is not controlled, it is possible that machine "runaway" might occur and the specimen would fail prematurely. For the case when the ratio of axial stress to tangential stress ($R = \sigma_A / \sigma_D$) is greater than or equal to unity, the axial strain is the control parameter and the computer program sends the command signal to the intensifier system to either increase or decrease the internal (or external) pressure so that a constant stress ratio will be maintained. On the other hand, when $R < 1$ the tangential strain will be the control parameter and the axial load is determined through the computer program. This procedure provides the stress ratio within 1% of a desired value.

III. DATA ANALYSIS

The biaxial test data presented here are in the form of effective stress-strain curves to illustrate strain hardening behavior and yield and failure stresses in the two-dimensional stress space. The data are presented in terms of the axial and tangential stresses σ_A and σ_T ; the corresponding strains are ϵ_A and ϵ_T .

Before discussing experimental results a definition of yield must be made. The definition used here is based on the square root of the second invariant of the stress deviator $\sqrt{J_2}$ plotted against the square root of the second invariant of the strain deviator $\sqrt{I_2}$, where in this study,

$$\sqrt{J_2} = \frac{1}{\sqrt{6}} [(\sigma_A - \sigma_T)^2 + (\sigma_T - \sigma_R)^2 + (\sigma_R - \sigma_A)^2]^{1/2}$$

$$\sqrt{I_2} = \frac{1}{\sqrt{6}} [(\epsilon_A - \epsilon_T)^2 + (\epsilon_T - \epsilon_R)^2 + (\epsilon_R - \epsilon_A)^2]^{1/2}$$

and σ_R is the radial stress and ϵ_R the radial strain.

The axial and tangential strains ϵ_A and ϵ_T are measured strain gage values on the outside surface of the tubular specimen. The radial strain ϵ_R through the wall of the tube is calculated using elasticity equations during elastic loading and the assumption of incompressible flow for the plastic components of strain after "yield."

The stresses are calculated from the axial load and the pressure. In the case of a tubular specimen subjected to axial load and internal pressure, the axial stress σ_A is

$$\sigma_A = \frac{F}{\pi(r_o^2 - r_i^2)}$$

where F is the axial load and r_o , r_i are the outer and inner radii and $r_o \geq r \geq r_i$ and the tangential stress σ_T is

$$\sigma_T = \frac{P_i r_i^2}{r_o^2 - r_i^2} \left(1 + \frac{r_o^2}{r^2} \right)$$

where P_i is the internal pressure.

In the case of a specimen subjected to axial load and external pressure, the vertical component of the hydraulic pressure acting on the surface of the fillet at both ends must be taken into consideration when the axial stress is calculated

$$\sigma_A = \frac{F}{\pi(r_o^2 - r_i^2)} + \frac{P_o (r_s^2 - r_o^2)}{(r_o^2 - r_i^2)}$$

where r_s = shoulder radius of the specimen (Figure 2) and P_o = external pressure. For specimen used in this study $r_s = 0.624$ inch.

The tangential stress

$$\sigma_T = - \frac{p_o r_o^2}{r_o^2 - r_i^2} \left(1 + \frac{r_i^2}{r^2} \right).$$

Once the deviatoric stress-strain curves are obtained, yielding is defined as the intercept of a line drawn parallel to the initial linear portion of the curve at 0.2% strain offset.

As it was discussed in the last section, all tests presented here are performed under proportional load conditions. In these tests, one strain rate was controlled to be constant while the second servo was used to maintain a constant stress ratio. Selection as to which direction would be maintained at constant strain rate was usually determined by the dominating stress. It will be shown below that tests performed under these control conditions also have a constant effective (or deviatoric) strain rate. The effective strain rate is defined as

$$\dot{\epsilon}_{eff}^P = \frac{1}{\sqrt{6}} \left[\left(\dot{\epsilon}_A^P - \dot{\epsilon}_T^P \right)^2 + \left(\dot{\epsilon}_R^P - \dot{\epsilon}_T^P \right)^2 + \left(\dot{\epsilon}_R^P - \dot{\epsilon}_A^P \right)^2 \right]^{1/2}$$

where the superscript P denotes plastic strain, and the dot indicates rate of change.

We further assume the Prandtl-Reuss flow rule

$$d\epsilon_A^P / S_A = d\epsilon_T^P / S_T = d\epsilon_R^P / S_R$$

or

$$\dot{\epsilon}_A^P / S_A = \dot{\epsilon}_T^P / S_T = \dot{\epsilon}_R^P / S_R$$

where S's are deviatoric stress components, and given as follows:

$$S_A = \sigma_A - 1/3 (\sigma_A + \sigma_T + \sigma_R) = 1/3 (2\sigma_A - \sigma_T)$$

$$S_T = \sigma_T - 1/3 (\sigma_A + \sigma_T + \sigma_R) = 1/3 (-\sigma_A + 2\sigma_T)$$

$$S_R = -1/3 (\sigma_A + \sigma_T)$$

with $\sigma_R = 0$.

The flow rule then becomes

$$\frac{\dot{\epsilon}_A^P}{1/3(2\sigma_A - \sigma_T)} = \frac{\dot{\epsilon}_T^P}{1/3(-\sigma_A + 2\sigma_T)} = \frac{\dot{\epsilon}_R^P}{-1/3(\sigma_A + \sigma_T)}.$$

If we define the stress ratio, $\beta \equiv \sigma_T/\sigma_A$, we have

$$\frac{\dot{\epsilon}_T^P}{\dot{\epsilon}_A^P} = \frac{2\sigma_T^{-\alpha} A}{2\sigma_A^{-\alpha} T} = \frac{2\beta-1}{2-\beta} \equiv \alpha$$

which gives $\dot{\epsilon}_T^P = \alpha \dot{\epsilon}_A^P$.

The assumption of incompressibility provides

$$\dot{\epsilon}_A^P + \dot{\epsilon}_T^P + \dot{\epsilon}_R^P = 0$$

and

$$\dot{\epsilon}_R^P = -(1+\alpha) \dot{\epsilon}_A^P.$$

Then substitution into the definition of effective strain rate yields

$$\dot{\epsilon}_{eff}^P = \dot{\epsilon}_A^P (1+\alpha+\alpha^2)^{1/2}.$$

Thus, for a selected proportional load path $\beta = \sigma_T/\sigma_A$ and a prescribed constant strain rate $\dot{\epsilon}_A^P$, the effective strain rate can be determined.

Two computer programs were written for the purpose of testing tubular specimens under biaxial stress states. The program "BIAX" was written to automatically record the axial load from the load cell, the pressure from the pressure transducer, and the axial and tangential strain values from the two sets of strain gages on the specimen. These values were stored on a magnetic tape through analog-to-digital converters. The program "BIANDV" was written to calculate the deviatoric stress, strain, and strain rate from the data recorded on the magnetic tape, and the results were printed on the line printer. A plotting routine was also included in "BIANDV," so the deviatoric stress-strain curve was automatically plotted on an X-Y recorder for each test.

The listing of these two programs are given in the Appendix.

IV. EXPERIMENTAL RESULTS AND DISCUSSIONS

The described biaxial testing technique was used to determine the yield surface of aluminum alloy 2014-T651 on the axial and tangential stress (σ_A, σ_T) plane. The tubular specimens were machined from 2- and 3-inch-thick rolled plates of 2014-T651 aluminum. Yield surface is shown in Figure 7. In the first quadrant (i.e., $\sigma_T > 0$, $\sigma_A > 0$) some specimens failed before reaching the 0.2% offset, but in the plastic region; in these cases, yield was defined as the maximum stresses reached during the test. The experimental data were also fitted to a quadratic equation with the least-square method

$$A\sigma_A^2 + B\sigma_A\sigma_T + C\sigma_T^2 + D\sigma_A + E\sigma_T = k^2$$

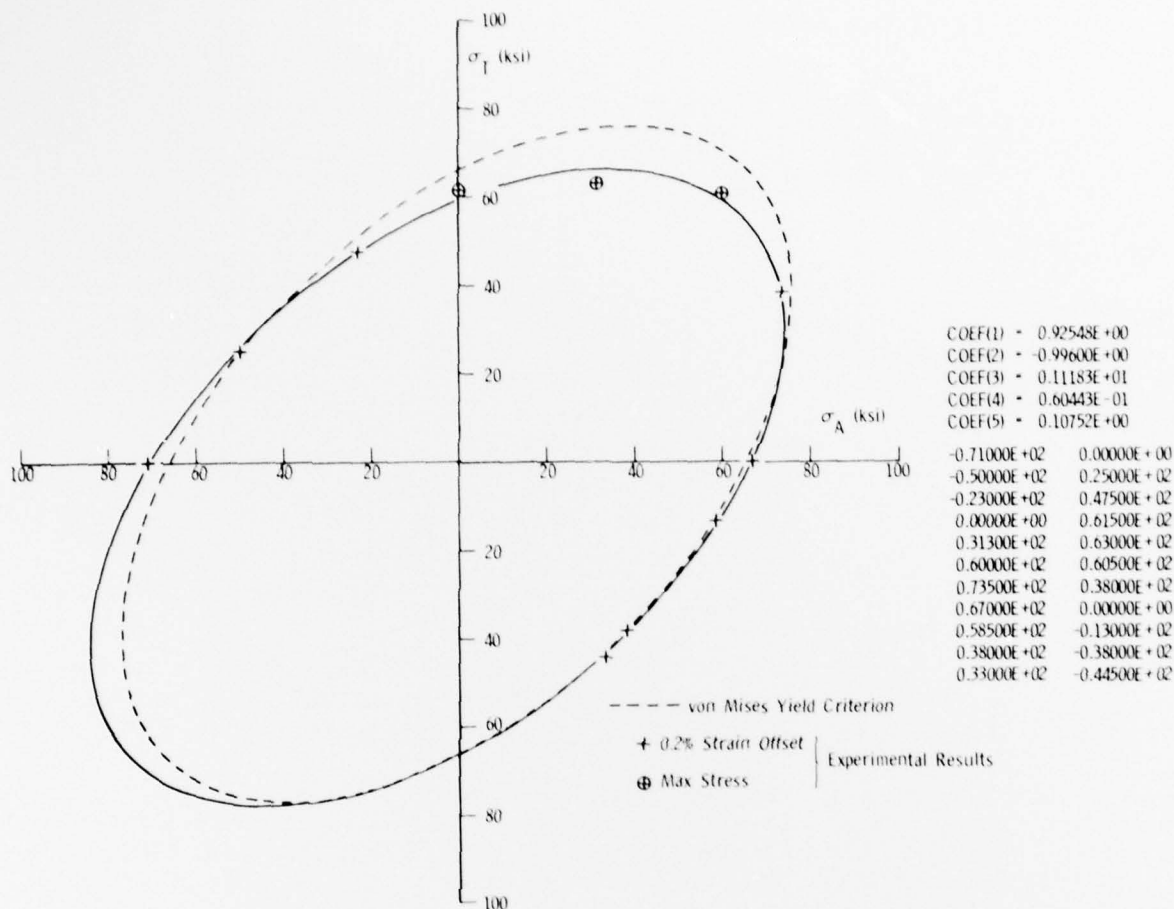


Figure 7. Yield surface for 2014-T651 aluminum alloy.

where k is the yield stress, 67 ksi, in simple tension. The coefficients are given as follows:

$$\begin{aligned} A &= 0.925 \\ B &= -0.996 \\ C &= 1.118 \\ D &= 0.060 \\ E &= 0.107 \end{aligned}$$

which indicates that the behavior of the material agrees reasonably well with the von Mises yield criteria. Furthermore, the experimental results presented in this report were compared with results obtained by other investigators, e.g., Reference 1, and they agree very well. This comparison serves as a verification of the testing techniques and data analysis procedure presented here.

1. REID, R. J., JONES, A. H., and GREEN, S. J. *Characterization of 2014-T651 Aluminum Alloy*. Terra Tek, Salt Lake City, Utah, Contract DAAG46-74-C-0019, Final Report, AMMRC CTR 74-68, November 1974.

APPENDIX. BIAX AND BIANDV PROGRAMS

05/8 FORTNAN IV 3.03

```

C      PROGRAM BIAX      TAKES IN 6 CH OF DATA
C
C      RATIO LESS THAN 1.0
C
C      LOAD CONTROL ON AXIAL DIRECTION
C
C      STRAIN CONTROL ON THETA DIRECTION
C
C      FIGURES LOAD COMMAND FROM PRESSURE READING
C
C      RATIO GREATER THAN 1.0 OR EQUAL
C
C      STRAIN CONTROL ON AXIAL DIRECTION
C
C      PRESSURE CONTROL ON THETA DIRECTION
C
C      FIGURES PRESSURE COMMAND FROM LOAD READING
C
0002      DIMENSION A(256),B(256),C(256),D(256),E(256),F(256),PLTBUF(256)
0003      DIMENSION DATHUF(50),STC(2)
C
C      SET UP DATA FILE
C
0004      DEFINE FILE 1(6,256,0,0)
0005      T=100
0006      Y=100
0007      DAI=-1
0010      CALL DTDA(1,DAI)
0011      DAI=-1
0012      CALL DTDA(2,DAI)
0013      WRITE(4,1059)
0014      1059 FORMAT(' IF EXTERNAL PRESSURE PUT ACTIVE GAGE IN
           1' ,/, ' DUMMY LOCATION = POSITIVE FEEDBACK ',/,)
0015      1000 FORMAT(' TEST NUMBER IS = ',S)
0016      1001 FORMAT(' DATE IS ',I3,I3,I3,/,)
0017      1002 FORMAT(' LOAD CAL = LBS = ',S)
0020      1003 FORMAT(' PRESS. CAL = PSI = ',S)
0021      1004 FORMAT(' INSIDE DIA. = ',S)
0022      1005 FORMAT(' OUTSIDE DIA. = ',S)
0023      1006 FORMAT(' TOTAL TEST TIME (SEC) ',F12.2)
0024      1007 FORMAT(' COMMAND SATURATED : TEST TERMINATED ')
0025      1008 FORMAT(' THETA STRAIN RATE CONSTANT ',/,
           1' AXIAL LOAD CONTROLLED BY AXIAL STRAIN FEEDBACK ',/)
0026      1009 FORMAT(' AXIAL LOAD : TYPE T OR C = ',S)
0027      1010 FORMAT(' STRESS RATIO = ',S)
0030      1011 FORMAT(' AXIAL LOAD SATURATED: TEST TERMINATED ')
0031      1012 FORMAT(' SAVING DATA ON UNIT 1 ')
0032      1020 FORMAT(' DO SET UP NOW : WHEN SET PUT GAGED SPECIMEN TO RUN ',/,
           1' HIT RETURN KEY',/)
0033      1021 FORMAT(' TYPE RETURN TO GO ')
0034      2000 FORMAT(2A4)
0035      2001 FORMAT(110)
0036      2002 FORMAT(1A1)

```

```

C
C      TEST NUMBER
C
0037      WRITE(4,1000)
0040      READ(4,2000) TEST1,TEST2
0041      CALL DATE(J1,J2,J3)
0042      WRITE(4,1001) J1,J2,J3
C
C      LOAD CALIBRATION
C
0043      WRITE(4,1002)
0044      READ(4,2001) ALCAL
C
C      PRESSURE CALIBRATION
C
0045      WRITE(4,1003)
0046      READ(4,2001) PCAL
C
C      INSIDE DIAMETER
C
0047      WRITE(4,1004)
0050      READ(4,2001) DI
0051      RI=DI/2.
C
C      OUTSIDE DIAMETER
C
0052      WRITE(4,1005)
0053      READ(4,2001) DO
0054      RO=DO/2.
C
C      STRESS RATIO    AXIAL/THETA
C
0055      WRITE(4,1010)
0056      READ(4,2001) RATIO
0057      21 CONTINUE
C
C      AXIAL LOAD = TENSION OR COMPRESSION
C
0060      WRITE(4,1009)
0061      READ(4,2002) AXIAL
0062      WRITE(4,1000)
0063      1060 FORMAT(' INTERNAL PRESSURE ? Y OR N = ',S)
0064      READ(4,2002) INT
0065      IF(INT.EQ.Y) GO TO 22
0066      WRITE(4,1001)
0067      1001 FORMAT(' SHOULDER DIA. = ',S)
0070      READ(4,2001) DS
0071      RS=DS/2.
C
C      SHOULDER AREA
C
0072      AS=3.14159*(RS**2-RO**2)
0073      CONE=PCAL*AS/511.
C
C      STRAIN GAGE INFORMATION
C

```

```

0074      22 00 401 I=1,2
0075      WRITE(4,1073) I
0076      1073 FORMAT(' *',I2,' RES. = ',S)
0077      READ(4,2001) RG
0100      WRITE(4,1072) I
0101      1072 FORMAT(' *',I2,' G.F. = ',S)
0102      READ(4,2001) GF
0103      WRITE(4,1074) I
0104      1074 FORMAT(' *',I2,' SHUNT = ',S)
0105      READ(4,2001) RS
0106      STC(I)=(1+RS/(RS+RG))/GF
0107      401 CONTINUE
C
C      COMPUTATIONS FOR EFFECTIVE STRAIN RATE
C      AND TEST TIME
C
0110      IF(RATIO.GE.1) BET=1./RATIO
0111      IF(RATIO.LT.1) BET=RATIO
0112      IF(AXIAL.NE.T.AND.INT.EQ.Y) BET=BET
0113      IF(AXIAL.EQ.T.AND.INT.NE.Y) BET=BET
0114      WRITE(4,1070)
0115      1070 FORMAT(' EFFECTIVE STRAIN RATE = ',S)
0116      READ(4,2001) EFF
0117      ALP=(2*BET-1)/(2.-BET)
0118      EOOT=EFF/SGHT(1.+ALP+ALP**2)
0120      WRITE(4,1071) EOOT
0121      1071 FORMAT(' CONTROLLED S.W. = ',F12.5)
0122      IF(RATIO.GE.1) TIME=STC(1)/EOOT
0123      IF(RATIO.LT.1) TIME=STC(2)/EOOT
0124      WRITE(4,1006) TIME
0125      11 WRITE(4,1021)
C
C      HIT RETURN KEY TO START TEST
C
0127      READ(4,2000) CONT
0130      INDA1=1
0131      IF(AXIAL.NE.T) INDA1=-1
0132      IF(AXIAL.NE.T) CALL DTGA(1,V)
0133      R=RG
0134      AA=3.14159*(R0**2-R1**2)
0135      TT=(R1**2*(R0**2+R**2))/(R**2*(R0**2-R1**2))
0136      IF(INT.NE.Y) TT=(R0**2*(R1**2+R**2))/(R**2*(R0**2-R1**2))
0137      CONV=20000*XCAL/(AA*TT*PCAL*511*RATIO)
0140      N0A0=H
0141      NPTS=256
0142      CALL CLRPLT(200,PLTHCF)
0143      KK=TIME/20.*1
0144      N0A=2000
0145      NPTS=(KK*2000)*H
0146      CALL NFALTM(DATBUF,50,H,B,NPTS)
0147      CR=(KK*2000)/TIME
C
C      INITIALIZE DATA TO ZERO
C
0150      00 300 I=1,256

```



```

0151      A(I)=0.0
0152      H(I)=0.0
0153      C(I)=0.0
0154      D(I)=0.0
0155      E(I)=0.0
0156      F(I)=0.0
0157      300 CALL PLOT(1,.05,.5)
0158      CALL CLOCK(N,CH)
0160      DO 301 I = 1,256
0161      IF(1.EQ.256) NOAAU=7
0162      DO 302 J=1,NOAAU
0163      CALL SSN(S,ISNSS)
0164      IF(ISNSS.EQ.1) GO TO 170
0165      AAA=0.0
0166      CC=0.0
0167      DO 303 K=1,KK
0168      CALL SSN(K,ISNSK)
0169      IF(ISNSK.EQ.1) GO TO 304
0170      IF(K.NE.KK) GO TO 304
0171
0172      C
0173      C      HARP ON CONTROL GAGE
0174      C
0175      IF(NATIO.LT.1) GO TO 340
0176      DAI=DAI+INDAI
0177      CALL DTGA(1,DAI)
0178      GO TO 304
0179      380 DAI=DAI+1
0180      CALL DTGA(2,DAI)
0181
0182      C
0183      C      SAMPLELINE ROUTINE
0184      C
0185      304 A(I)=AOB(X)
0186      B(I)=AOB(X)
0187      C(I)=AOB(X)
0188      D(I)=AOB(X)
0189      E(I)=AOB(X)
0190      F(I)=AOB(X)
0191      IF(INT.EQ.Y) C(I)=-C(I)
0192      IF(INT.EQ.Y) D(I)=-D(I)
0193      IF(INT.NE.Y) A(I)=(A(I)+XLCAL/511.+ABS(C(I))*CONE)*511./XLCAL
0194      AAA=AAA+A(I)
0195      CC=CC+C(I)
0196      303 CONTINUE
0197
0198      C
0199      C      AVERAGE DATA BEFORE STORAGE
0200      C
0201      A(I)=AAA/KA
0202      C(I)=CC/KA
0203      IF(NATIO.GE.1.0) GO TO 383
0204
0205      C
0206      C      COMMAND TO MSHH
0207      C
0208      SIGZ=(XLCAL/(AA))
0209      SIGT=ABS(C(I)+PCAL*IT)*NATIO
0210      DAI=SIGT/SIGZ*4,

```

```

0224      IF (AXIAL.NE.T) OA1=-OA1
0225      IF (INT.NE.Y) OA1=OA1+C(1)*CONE*511./XLCAL*4.
0226      IF (ABS(OA1).GT.2047) GO TO 920
0227      CALL SSW(0,ISNS0)
0230      IF (ISNS0.EQ.1) GO TO 313
0231      CALL DTDA(1,OA1)
0232      GO TO 313

C
C      COMMAND TO INTENSIFIER
C
0233      303 OA2=ABS(A(1)*CONV)
0234      IF (OA2.GT.2047) GO TO 920
0235      CALL SSW(0,ISNS0)
0236      IF (ISNS0.EQ.1) GO TO 313
0237      CALL DTDA(2,OA2)
0240      313 CONTINUE
0241      302 CONTINUE

C
C      PLOT AXIAL VS THETA STRESS ON SCOPE
C
0242      CALL PLOT(1,A(1)*1.5/1024,+.65,C(1)/1024,+.5,1)
0243      301 CONTINUE
0244      GO TO 100
0245      800 CONTINUE
0246      WRITE(4,1011)
0247      GO TO 100
0250      900 WRITE(4,1017)
0251      100 CONTINUE

C
C      DUMP DATA TO TAPE
C
0252      WRITE(4,1010)
0253      L=1
0254      WRITE(1*L) A
0255      L=2
0256      WRITE(1*L) B
0257      L=3
0260      WRITE(1*L) C
0261      L=4
0262      WRITE(1*L) D
0263      L=5
0264      WRITE(1*L) E
0265      L=6
0266      WRITE(1*L) F
0267      WRITE(4,1008)
0270      1008 FORMAT(' RETURN TO END PROG. & REMOVE SOFTWARE CLAMP')
0271      READ(4,2001) GO
0272      CALL WOPEN(0)
0273      CALL WOPEN(2)
0274      103 CALL EXIT

```

05/6 FORTRAN IV 3.23

```

C PROGRAM HIANOV ANALYSIS OF A CH OF HIAK
C
C RATIO LESS THAN 1.
C
C STRAIN CONTROL ON AXIAL DIRECTION
C
C STRAIN CONTROL ON THETA DIRECTION
C
C FIGURES LOAD COMMAND FROM PRESSURE READING
C
C RATIO GREATER THAN 1.
C
C STRAIN CONTROL ON AXIAL DIRECTION
C
C PRESSURE CONTROL ON THETA DIRECTION
C
C FIGURES PRESSURE COMMAND FROM LOAD READING
C
C
C COMPUTES DEVIATORIC STRESS-STRAIN CURVE
C
C PLOTS DEVIATORIC STRESS-STRAIN CURVE
C
C
C USES SUBROUTINE XYREC
C
C
0002 DIMENSION A(256),B(256),C(256),D(256),E(256),F(256)
0003 DIMENSION NS(4),NG(4),STC(4)
0004 COMMON/CA,B,C,D,E,F
0005 DEFINE FILE 1(6,256,0,0)
0006 Y=1MY
C
C READ DATA FROM TAPE
C
0007 L=1
0010 READ(1*L) A
0011 L=2
0012 READ(1*L) H
0013 L=3
0014 READ(1*L) C
0015 L=4
0016 READ(1*L) D
0017 L=5
0020 READ(1*L) E
0021 L=6
0022 READ(1*L) F
C
C CROSS PLOT DATA ON SCOPE
C
C RETURN TO EXIT FROM PLT6
0023 CALL PLT6
0024 1000 FORMAT(' TEST NUMBER IS ',2A4)
0025 1001 FORMAT(' DATE IS ',13,13,13,13,13,13)
0026 1002 FORMAT(' LOAD CAL = LBS ')

```

THIS PAGE IS BEST QUALITY PRACTICABLE
FROM COPY FURNISHED TO DDC

```

0027 1003 FORMAT(' PRESS. CAL = PSI ')
0030 1004 FORMAT(' INSIDE DIA. ')
0031 1005 FORMAT(' OUTSIDE DIA. ')
0032 1006 FORMAT(' TOTAL TEST TIME (SEC) ')
0033 1010 FORMAT(' STRESS RATIO ')
0034 1011 FORMAT(' INTERNAL PRESSURE ? ',3)
0035 1021 FORMAT(' TYPE RETURN TO GO ')
0036 2000 FORMAT(240)
0037 2001 FORMAT(110)
0040 2002 FORMAT(141)

C
C INPUT CALIBRATION VALUES
C
0041 WRITE(4,1000)
0042 READ(4,2000) TEST1,TEST2
0043 CALL DATE(J1,J2,J3)
0044 WRITE(4,1001) J1,J2,J3
0045 WRITE(4,1002)
0046 READ(4,2001) PCAL
0047 WRITE(4,1003)
0048 READ(4,2001) PCAL
0049 WRITE(4,1004)
0050 READ(4,2001) DI
0051 READ(4,2001) DI
0052 WRITE(4,1005)
0053 READ(4,2001) DO
0054 READ(4,2001) DO
0055 WRITE(4,1006)
0056 READ(4,2001) DO
0057 READ(4,2001) DO
0058 WRITE(4,1007)
0059 READ(4,2002) INT
0060 WRITE(4,1008)
0061 READ(4,2001) RATIO
0062 WRITE(4,1009)

C
C INPUT STRAIN GAGE INFORMATION
C
0065 READ(4,2001) TIME
0066 WRITE(4,1010)
0067 READ(4,2001) GF
0070 1101 FORMAT(' GAGE FACTOR ')
0071 DO 12 I=1,4
0072 WRITE(4,1102) I
0073 1102 FORMAT(' *',I2,' GAGE RES. ')
0074 12 READ(4,2001) RG(I)
0075 DO 13 I=1,4
0076 WRITE(4,1103) I
0077 1103 FORMAT(' *',I2,' SHUNT RES. ')
0080 13 READ(4,2001) RS(I)
0081 DO 14 I=1,4
0082 14 STC(I)=(1-RS(I))/(RS(I)+RG(I))/GF*100

C
C SETUP FOR DETERMINING PRINCIPLE STRESSES
C
0103 A=(R1**2+PCAL/511)/(R0**2-R1**2)
0104 IF(INI.NE.V) A=(R0**2+PCAL/511)/(R0**2-R1**2)

```

```

0105      SIG1C=(ALCAL/511.)/((3.14159*(RO**2-NI**2))/1000.
0106      SIG2C=AAA*(1.+(RO**2/R**2))/1000.
0107      IF(INT,NE,Y) SIG2C=AAA*(1-NI**2/R**2)/1000.
0110      SIG3C=AAA*(1.-(RO**2/R**2))/1000.
0111      IF(INT,NE,Y) SIG3C=AAA*(1-(NI**2/R**2))/1000.
C
C      DETERMINE IF TWO AXIAL OR TWO TRANSVERSE
C      GAGES DIFFER BY MORE THAN 5% OF CALIBRATION VALUE
C
0112      DO 20 I=1,250
0113      E1=B(1)*STC(1)
0114      E2=B(1)*STC(2)
0115      E3=B(1)*STC(3)
0116      E4=B(1)*STC(4)
0117      DE1=(ABS(E1)-ABS(E3))
0120      DE2=(ABS(E2)-ABS(E4))
0121      DE1=.5*STC(1)*511
0122      DE2=.5*STC(2)*511
0123      IF(DE1,GT,DE1) GO TO 21
0124      IF(DE2,GT,DE2) GO TO 22
0125      20 CONTINUE
0126      GO TO 27
0127      21 WRITE(4,1104) I
0130      GO TO 27
0131      22 WRITE(4,1105) I
0132      1104 FORMAT(' AXIAL GAGES DIFFER THAN MORE THAN 5% AT PT #',13)
0133      27 WRITE(4,1141)
0134      1141 FORMAT(' MODULUS X10**6 PSI ')
0135      READ(4,2011) EMOU
0136      WRITE(4,1142)
0137      1142 FORMAT(' POISSONS RATIO = ')
0140      READ(4,2011) U
0141      EMOU=EMOU*1.000000.
0142      1105 FORMAT(' THETA GAGES DIFFER MORE THAN 5% AT PT #',13)
0143      23 J=1
0144      DO 30 I=1,250
0145      IF(I,NE,J,OR,ISNSC,EG,1) GO TO 40
0146      WRITE(3,5000)
0147      5000 FORMAT(1H1)
0150      WRITE(3,1000) TEST1,TEST2
0151      WRITE(3,1001) J1,J2,J3
0152      WRITE(3,1106)
0153      1106 FORMAT(' AXIAL THETA RADIAL AXIAL THETA',
1' TIME STRAIN DEV DEV DEV ')
0154      WRITE(3,1007)
0155      1007 FORMAT(' STRESS STRESS STRESS STRAIN STRAIN',
1' SEC RATE STRESS STRAIN RATE ',//)
0156      J=J+52
C
C      COMPUTE PRINCIPLE & DEVIATIONIC STRESSES
C
0157      40 S1=SIG1C*A(1)
0158      S2=SIG2C*C(1)
0159      S3=SIG3C*C(1)
0160      E1=((B(1)*STC(1))+(E(1)*STC(3)))/2./511.

```

THIS PAGE TO BEST QUALITY PRACTICABLE
FROM COPY FURNISHED TO DDC


```

0163      E2=((O(1)*STC(2))*(F(1)*STC(4)))/2./511.
0164      T1=(S2-S3)**2
0165      T2=(S3-S1)**2
0166      T3=(S1-S2)**2
0167      SIGO=(1/SGHT(6))*SGHT(T1+T2+T3)

C
C      COMPUTE RADIAL STRAIN AND DEVIATORIC STRAIN
C
0170      E(I)*SIGO
0171      E3=1/EMOD*(S3*1000.-O*(S1+S2)*1000.)*100.
0172      F1=((E2-E3)/100.)*2
0173      F2=((E3-E1)/100.)*2
0174      F3=((E1-E2)/100.)*2
0175      EEO=(1/SGHT(6))*SGHT(F1+F2+F3)*100.
0176      F(I)*EEO
0177      TM=(TIME/256)*1
0178      T1=TIME/256
0200      CALL SS*(O,ISNS,O)
0201      IF(1,GT,15) GO TO 50
0202      SH=0.0
0203      WRITE(3,1130) S1,S2,S3,E1,E2,TM,SH,SIGO,EEO
0204      GO TO 32
0205
50      IF(NATIO,LT,1) SH*(O(I)-O(I-15))*STC(2)/(T1*15)/100./511.
0206      IF(NATIO,GE,1) SH*(B(1)-B(I-15))*STC(1)/(T1*15)/100./511.
0207      IF(NATIO,LT,1) GO TO 37
0210      1130 FORMAT(3F8.2,2F8.3,F8.2,F4.5,2F8.2,F4.5)
0211      37 CONTINUE
0212      IF(S2,EG,0.0) GO TO 888
0213

C
C      DETERMINE EFFECTIVE STRAIN RATE
C
0214      RAT=S1/S2
0215      888 IF(S2,EG,0.0) RAT=S1/99999
0216      IF(NATIO,GE,1) BET=1./NATIO
0217      IF(NATIO,LT,1) BET=NATIO
0218      IF(RAT,LT,0) BET=-BET
0221      ALP=(2.*BET-1.)/(2.*BET)
0222      DEV3=SH*SGHT(1.+ALP*ALP**2)
0223      IF(ISNS,EG,1) GO TO 32
0224      WRITE(3,1130) S1,S2,S3,E1,E2,TM,SH,SIGO,EEO,DEV3
0225      32 CONTINUE
0226      32 WRITE(3,5000)

C
C      PLOT DEVIATORIC STRESS-STRAIN CURVE ON
C      X-Y RECORDER
C
0227      WRITE(4,1143)
0230      1143 FORMAT(' X-Y PLOT X=CH#0, Y=CH#1, YES OR NO ')
0231      READ(4,2002) PLOT
0232      Y=1-PY
0233      IF(PLOT,NE,Y) CALL EXIT
0234      CALL XYREC
0235      STOP
0236      END

```

DISTRIBUTION LIST

No. of Copies	To
	Office of Deputy Under Secretary of Defense for Research and Engineering, The Pentagon, Washington, D.C. 20301
1	ATTN: Mr. J. Persh, Staff Specialist for Materials and Structures, Room 3D1089
12	Commander, Defense Documentation Center, Cameron Station, Building 5, 5010 Duke Street, Alexandria, Virginia 22314
	Commander, U.S. Army Materiel Development and Readiness Command, 5001 Eisenhower Avenue, Alexandria, Virginia 22333
1	ATTN: DRCLDC, Mr. R. Zentner
	Director, Ballistic Missile Defense Advanced Technology Center, P.O. Box 1500, Huntsville, Alabama 35807
1	ATTN: ATC-M, Dr. D. Harmon
1	ATC-M, Mr. M. Whitfield
1	ATC-M, Mr. J. Papadopoulos
	Commander, Ballistic Missile Defense Systems Command, P.O. Box 1500, Huntsville, Alabama 35807
1	ATTN: BMDSC-TEN, Mr. N. J. Hurst
	Director, Defense Nuclear Agency, Washington, D.C. 20305
1	ATTN: SPAS, Mr. J. F. Moulton, Jr.
	Deputy Chief of Staff, Research, Development and Acquisition, Department of the Army, Washington, D.C. 20310
1	ATTN: DAMA-CSS, Dr. J. Bryant
	Director, U.S. Army Ballistic Research Laboratory, Aberdeen Proving Ground, Maryland 21005
1	ATTN: Mr. J. Meszaros
1	Dr. G. L. Moss
	Commander, U.S. Army Missile Research and Development Command, Huntsville, Alabama 35809
1	ATTN: Dr. R. Rhodes
1	Dr. S. Smith
	Commander, Harry Diamond Laboratories, 2800 Powder Mill Road, Adelphi, Maryland 20783
1	ATTN: DELHD-RBF, Dr. R. Oswald
1	DELHD-NP, Dr. F. Wimenitz
	Commander, Picatinny Arsenal, Dover, New Jersey 07801
1	ATTN: Mr. M. Allen
1	Mr. M. Weinstein
1	Mr. B. Frank

No. of Copies	To
1	Commander, U.S. Army Combat Development Command, Institute of Nuclear Studies, Fort Bliss, Texas 79916 ATTN: Technical Library
1	Commander, Air Force Materials Laboratory, Air Force Systems Command, Wright-Patterson Air Force Base, Ohio 45433 ATTN: LNE/1LT T. Curci
1	LNC/Dr. D. Schmidt
1	Department of the Navy, Naval Ordnance Systems Command, Washington, D.C. 20360 ATTN: ORD-03331, Mr. M. Kinna
1	Commander, Naval Surface Weapons Center, Silver Spring, Maryland 20910 ATTN: Mr. L. Gowen
1	Mr. F. Koubek
1	Los Alamos Scientific Laboratory, P.O. Box 1663, Los Alamos, New Mexico 87544 ATTN: GMX-6, Dr. J. W. Taylor
1	Space and Missile Systems Organization, P.O. Box 92960, World Way Postal Center, Los Angeles, California 90009 ATTN: RSSE/LTC J. McCormack
1	RSSE/CAPT J. Bolen
1	Sandia Laboratories, P.O. Box 5800, Albuquerque, New Mexico 87115 ATTN: Dr. D. Munson
1	Dr. W. Herrmann
1	Dr. L. D. Bertholf
1	Dr. B. Butcher
1	Dr. J. Lipkin
1	Aerospace Corporation, P.O. Box 92957, Los Angeles, California 90009 ATTN: Dr. R. Cooper
1	Dr. W. Barry
1	AVCO Corporation, Government Products Group, 201 Lowell Street, Wilmington, Massachusetts 01997 ATTN: Dr. W. Reinecke
1	Mr. P. Rolincik
1	Effects Technology, Inc., P.O. Box 30400, Santa Barbara, California 93105 ATTN: Dr. R. Wengler
1	Mr. J. Green
1	Mr. E. Steele
1	Fiber Materials, Inc., Biddeford Industrial Park, Biddeford, Maine 04005 ATTN: Mr. Maurice Subilia, Jr.
1	Mr. L. Landers
1	Mr. G. Williams

No. of
Copies

To

General Electric Company, Valley Forge Space Technology Center,
P.O. Box 8555, Philadelphia, Pennsylvania 19101
1 ATTN: Mr. K. Hall
1 Mr. R. Sullivan
1 Mr. J. Brazel

Kaman Sciences Corporation, P.O. Box 7463, Colorado Springs, Colorado 80933
1 ATTN: Mr. F. Shelton

Ktech, 911 Pennsylvania Avenue, N.E., Albuquerque, New Mexico 87110
1 ATTN: Dr. D. Keller

Lockheed Missiles and Space Company, P.O. Box 504, Sunnyvale,
California 94088
1 ATTN: Mr. D. Aspinwall

Martin Marietta Aerospace, P.O. Box 5837, Orlando, Florida 32805
1 ATTN: Mr. J. Potts
1 Mr. L. Kinnaird
1 Mr. F. Koo

McDonnell Douglas Corporation, 5301 Bolsa Avenue, Huntington Beach,
California 92647
1 ATTN: Mr. N. Campbell

Prototype Development Associates, Inc., 1740 Garry Avenue, Suite 201,
Santa Ana, California 92705
1 ATTN: Dr. McDonald
1 Mr. J. Schutzler

R&D Associates, 525 Wilshire Blvd., Santa Monica, California 90025
1 ATTN: Dr. A. Field

Southwest Research Institute, 8500 Culebra Road., San Antonio, Texas 78206
1 ATTN: Mr. A. Wenzel

Stanford Research Institute, 333 Ravenswood Avenue, Menlo Park,
California 90250
1 ATTN: Dr. D. Curran
1 Dr. L. Seaman

Stanford University, Department of Applied Mechanics, Stanford,
California 94305
1 ATTN: Prof. E. H. Lee

Director, Army Materials and Mechanics Research Center,
Watertown, Massachusetts 02172
2 ATTN: DRXMR-PL
1 DRXMR-AG-MD
3 Authors

<p>Army Materials and Mechanics Research Center, Watertown, Massachusetts 02172 BIAxIAL TESTING TECHNIQUES OF THIN-WALLED TUBULAR SPECIMENS - James H. Rainey, Ronald A. Swanson, and Shun-Chin Chou</p> <p>Technical Report AMMRC TR 78-37, September 1978, 22 pp., illus., D/A Project 8X363304D215, AMMCS Code 633304.21500.03</p> <p>A biaxial testing technique of a thin-walled tubular specimen is described. The specimen was loaded with a combination of axial tension/compression and internal/external pressure. The technique can be used to test tubular specimens of any material. The yield surface of 2014-T651 aluminum was determined to illustrate the testing technique and data analysis.</p>	<p>AD</p> <p>UNCLASSIFIED UNLIMITED DISTRIBUTION</p> <p>Key Words Biaxial stresses Mechanical tests Yield</p>
---	---

<p>Army Materials and Mechanics Research Center, Watertown, Massachusetts 02172 BIAxIAL TESTING TECHNIQUES OF THIN-WALLED TUBULAR SPECIMENS - James H. Rainey, Ronald A. Swanson, and Shun-Chin Chou</p> <p>Technical Report AMMRC TR 78-37, September 1978, 22 pp., illus., D/A Project 8X363304D215, AMMCS Code 633304.21500.03</p> <p>A biaxial testing technique of a thin-walled tubular specimen is described. The specimen was loaded with a combination of axial tension/compression and internal/external pressure. The technique can be used to test tubular specimens of any material. The yield surface of 2014-T651 aluminum was determined to illustrate the testing technique and data analysis.</p>	<p>AD</p> <p>UNCLASSIFIED UNLIMITED DISTRIBUTION</p> <p>Key Words Biaxial stresses Mechanical tests Yield</p>
---	---

<p>Army Materials and Mechanics Research Center, Watertown, Massachusetts 02172 BIAXIAL TESTING TECHNIQUES OF THIN-WALLED TUBULAR SPECIMENS - James H. Rainey, Ronald A. Swanson, and Shun-Chin Chou</p> <p>Technical Report AMMRC TR 78-37, September 1978, 22 pp., illus., D/A Project 8X363304D215, AMCMS Code 633304.21500.03</p> <p>A biaxial testing technique of a thin-walled tubular specimen is described. The specimen was loaded with a combination of axial tension/compression and internal/external pressure. The technique can be used to test tubular specimens of any material. The yield surface of 2014-T651 aluminum was determined to illustrate the testing technique and data analysis.</p>	<p>AD</p> <p>UNCLASSIFIED UNLIMITED DISTRIBUTION</p> <p>Key Words</p> <p>Biaxial stresses Mechanical tests Yield</p>
---	--

<p>Army Materials and Mechanics Research Center, Watertown, Massachusetts 02172 BIAXIAL TESTING TECHNIQUES OF THIN-WALLED TUBULAR SPECIMENS - James H. Rainey, Ronald A. Swanson, and Shun-Chin Chou</p> <p>Technical Report AMMRC TR 78-37, September 1978, 22 pp., illus., D/A Project 8X363304D215, AMCMS Code 633304.21500.03</p> <p>A biaxial testing technique of a thin-walled tubular specimen is described. The specimen was loaded with a combination of axial tension/compression and internal/external pressure. The technique can be used to test tubular specimens of any material. The yield surface of 2014-T651 aluminum was determined to illustrate the testing technique and data analysis.</p>	<p>AD</p> <p>UNCLASSIFIED UNLIMITED DISTRIBUTION</p> <p>Key Words</p> <p>Biaxial stresses Mechanical tests Yield</p>
---	--